

This model features electric operation of height and tilt facility with manually adjustable backrest and individually adjustable leg rests to aid positioning of the patient for podiatry procedures. The backrest has auto levelling armrests to provide support and comfort to the patient in sitting positions. The backrest and leg sections can also be adjusted to fully support the patient in prone or supine positions for biomechanics or other procedures. A wide range of options are available to personalise the couch for various procedures and to meet individual user requirements.


Model 93CDT Electric
Colour - Sapphire Blue



Standard Features

- UK designed & manufactured
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Safe working load & lift capacity 250kg
- Electric height & tilt operation
- Manual adjustment of backrest angle
- Manual adjustment of each leg section
- Extendable heel pads
- Self-levelling armrests
- Hardwearing antimicrobial vinyl
- Comfort foam
- Retractable wheel base
- Adjustable foot for uneven floors
- Hoist clearance 13cm (5")

Upholstery Dimensions

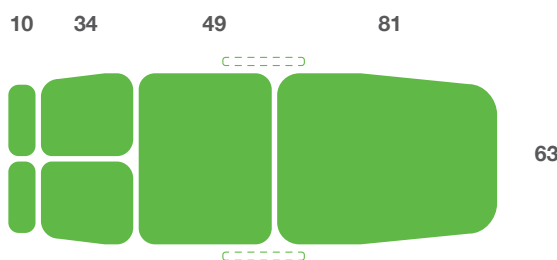

Overall Dimensions

Length: 176cm (69") approx
Width: 63cm (25") approx
Height: 47 - 97cm approx
(18 - 38") approx

Model 93CDT Dual Electric

Options & Accessories



Breathing
Hole



Paper Roll
Holder



Footswitch



Motorised
Backrest



D Cushion



Operator's
Chair

Operational Parameters

Backrest range +85° to -15°



Trendelenburg Tilt -20°

Foot section range +10° to -90°



Model 93CDT Dual Electric.

Upholstery Vinyl - Phthalate Free

- Flame retardant to CRIB 5
- Independantly tested to verify effectiveness at reducing the spread of MRSA, E.coli and other similar bacteria and fungi
- Rated at 1,000,000 Martindale rubs without damaging or marking the top layer, making it suitable for all high traffic, heavy duty areas



Jasmine



Battleship



Colbalt



Cocoa



Grape



Mulled Wine



Pillarbox



Gingersnap



Candy



Applemint



Sapphire



Atlantic Blue



Lupin



Denim



Cool Blue



Rainforest



Wasabi



Citrus Green



Marigold



Almond

Product Benefits

- UK designed & manufactured with long-term support
- Lifetime frame and 5 year parts guarantee providing longevity of product life and value for money
- Safe working load of 250kgs to cater for very heavy patients
- Simple push button electric height and tilt operation
- Trendelenburg tilt improves patient positioning in sitting as well as providing a recovery position
- Manual adjustment of backrest angle to support the patient in sitting or lying positions.
- Auto levelling armrests for patient support & comfort in sitting or semi-sitting positions.
- The armrests fold back to enable side access
- The leg sections are manually adjustable from 10° to -90° to aid positioning of the lower limbs
- Retractable wheel base system for ease of transport and positioning
- Adjustable foot for stability on uneven floor surfaces
- Option of footswitch control for height & tilt operation
- Easy to clean antimicrobial upholstery vinyl
- Hoist clearance below frame 13cm (5")

Technical Standards

- Safe Working Load & Lifting Capacity tested in accordance with the rating on all models
- MD Classification (Class 1) EU MDR 2017/745 Annex VIII Class 1 / UK MDR 2002 (SI 2002 No 618, as amended)
- Electrical Standard - IEC 60601-1
- IP Rating IPX6/IPX4 depending on options selected
- Vinyl flame retardant to CRIB 5 (BS 5852 / BS7176)
- Foam flame retardant to BS 5852: Part 2: 1982, Ignition Source 5 (CRIB 5)
- UKCA and CE Marked Medical Device
- Plinth Medical Ltd is externally audited by LRQA to ISO 9001:2015

Plinth Medical Ltd.
Barric Lane
Occold
Eye
Suffolk, IP23 7PX

t: +44 (0) 1449 767887
e: sales@plinthmedical.com

www.plinthmedical.com

Plinth | Medical



Version 1 January 2024 E&OE