



FACE SHIELD

PRODUCT CODE: TEC-FACEVISOR _ 2

DISPOSABLE FACE SHIELD


WITH SOFT ELASTICATED STRAP

LIGHTWEIGHT FACE SHIELD

Tecman's lightweight face shield with anti-fog properties protects the face and eyes from splashing or spraying of infectious or hazardous substances. The disposable face shields are manufactured from optically clear polyester for maximum visibility, whilst the soft elastic strap and PU foam headband provide additional, long-term comfort and durability.

The product can be worn with or without additional face protection, such as prescription or protective eyewear and face masks.

Our registered design face shield is manufactured using a patent-protected process.

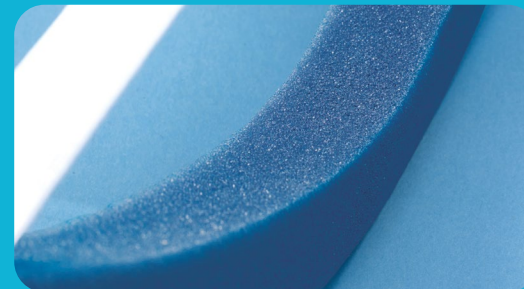

KEY FEATURES

- ❖ Anti-fog properties*
- ❖ Optically clear visor for maximum visibility
- ❖ PU foam headband for long-term comfort
- ❖ Soft elasticated strap
- ❖ Lightweight design
- ❖ Protects the face and eyes from splashing or spraying
- ❖ Can be worn with or without eyewear protection and face masks
- ❖ EN 166:2001 category III certificated (EU) 2016/425
- ❖ Community Design No. 008192637

SPECIFICATIONS

SPEC	DETAIL
PRODUCT COLOUR	OPTICALLY CLEAR (NO TINT)
PRODUCT TYPE	FACE SHIELD
SIZE	330MM X 220MM
OPTICAL CLASS	1 (+/-0.06)
FIELD OF USE	DROPLETS OR SPLASH OF LIQUIDS
APPLICATION	HEALTHCARE, PPE (NOT FOR INDUSTRIAL USE)
STANDARDS	EN166:2001 CERT REF: 2777/14610-02/E00-00 - CATEGORY III - COMPLEX DESIGN cl 6.1-6.3, 7.1.1, 7.1.2.1-7.1.2.3, 7.1.3, 7.1.5.1, 7.1.7 & 7.2.4
SHELF LIFE	2 YEARS
STORAGE CONDITIONS	STORE IN ORIGINAL PACKAGING OUT OF DIRECT SUNLIGHT AT 23°C AND 50% RH

MATERIAL DATA

MATERIAL	DESCRIPTION	SIZE
VISOR FILM	OPTICAL GRADE PET	0.18MM THICK
FOAM HEADBAND	SOFT GRADE PU	25MM THICK
ELASTIC STRAP	MEDICAL GRADE, LATEX FREE	19MM WIDE

INSTRUCTIONS FOR USE


Carefully remove the protection film from the rear of the visor before use. To put on and take off products ensure the elastic is wrapped securely around the head with foam against the forehead. Only wear products of a suitable size that are correctly adjusted. Products which are loose may not stay in place and if too tight will be uncomfortable. The size of this product is available in one standard size with an elastic headband.

ADDITIONAL SPECIFICATIONS

Latex Free / Medical grade materials / Single use, cleaning not recommended / Can be used for prolonged periods as it does not create distortion or fatigue / Visor suitable for use with prescription lenses and surgical masks / Always check before using / Manufactured in the UK

*Not compliant with EN166: 2001 7.3.2 Resistance to Fogging

DOC252, ISSUE 2, FEB 2021

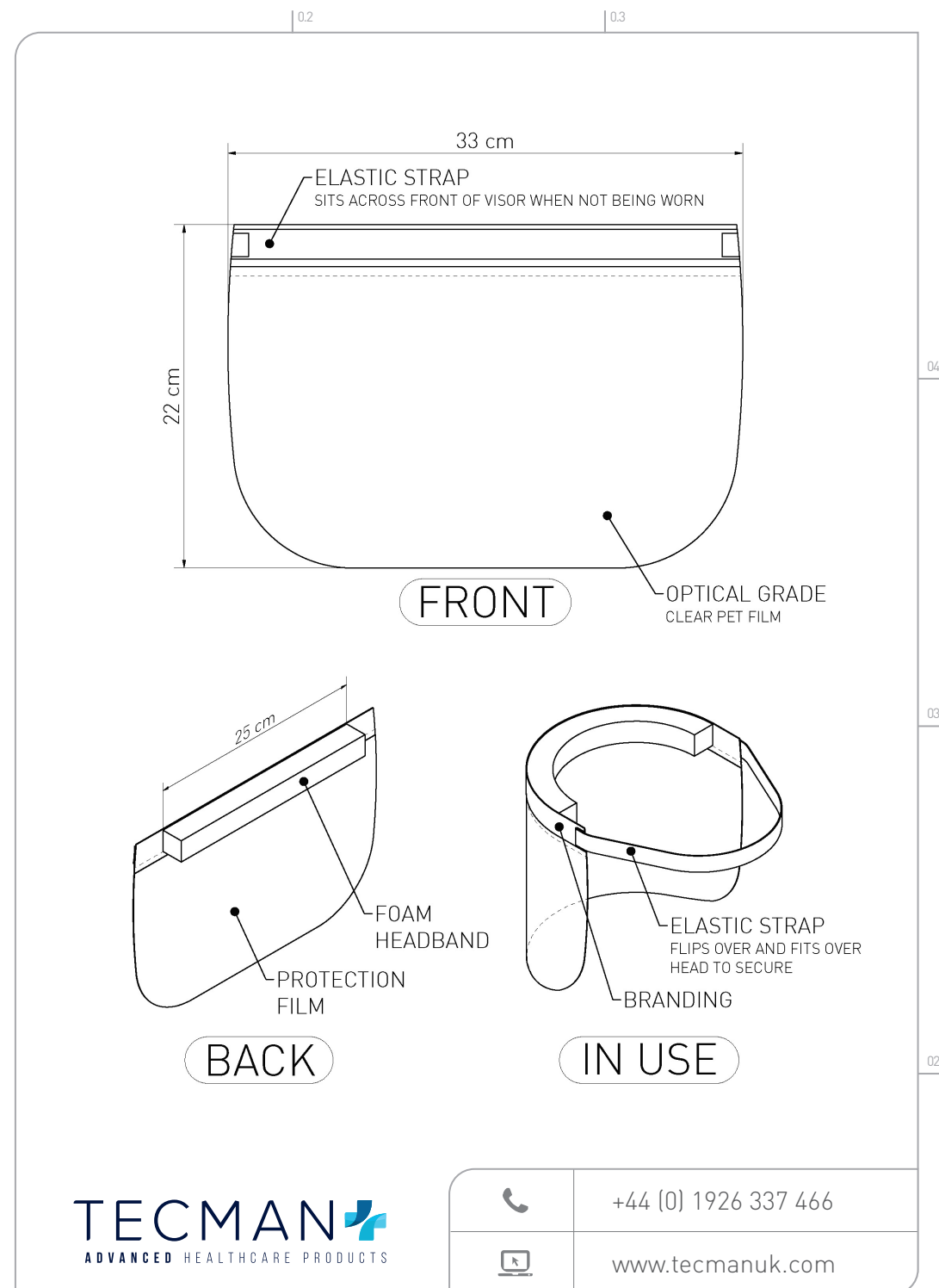

TECMAN
ADVANCED HEALTHCARE PRODUCTS



+44 (0) 1926 337 466



www.tecmanuk.com